国家高新技术企业 ISO 13485 YY/T 0287











高品质复合材料机箱 Melamine unit box

整个机箱采用美耐皿树脂, 光滑亮丽, 不易变色, 表面 易于清洗消毒。

The unit box is made of melamine resin, which is glossy and bright, not easy to discolor, easy to clean and disinfect.

无贴纸一体副控 No-sticker integrated assistant control

- 1.防止液体从贴纸的缝隙进入线路板,提高电路可靠性,侧边带拉手便于牙医操作。 2.强弱吸可延时5秒关闭,以确保强弱吸管内污物被完全排空。
- 3.强弱吸手柄采用铝合金和硅胶组合构成,可拆卸清洗消毒,重量轻,密封性好,可通 过手柄流量开关,对流量进行调节。
- 1. Prevent the liquid from entering the circuit board through the gap of the sticker, improve the reliability of the circuit. Side handle makes operation easier.
- 2. The strong suction and saliva ejector can be turned off with a delay of 5 seconds to ensure that the dirt in tubes is completely drained.
- 3. The strong suction and saliva ejector handles are made of a combination of aluminum alloy and silicone, light weight and good sealing, can be removed for cleaning and disinfection. The flow can be adjusted by the handle flow switch.





Trolley-type doctor table (for option)

- 1.推车式器械盘改善牙椅承载情况,使整机更稳定,消除器械臂给患者带来的压迫感。
- 2.底部采用铝合金压铸工艺,提高结构稳定性,不变形、不生锈。
- 3.倾斜式的按键面板设计,更具人性化避免牙医的误操作。
- 4.观片灯位置可升级至彩色显示屏,控制更轻松,显示更清晰(选配)。
- 1. Trolley-type doctor table reduces the burden of the whole unit and makes the chair more stable, eliminates the oppressive feeling from patients caused by instrument arm.
- 2. The base adopts aluminum alloy die-casting technology to improve the structural stability without deformation and rust.
- 3. Tiling design of key board is more user-friendly and good for accurate operation.
- 4. In stead of X-film viewer, the screen can be upgraded into color display for easier control and clearer vision (for option).

创新引领未来 Innovation leads the future



旋转陶瓷痰盂

Rotary ceramic cuspidor

可旋转的陶瓷痰盂便于病人吐痰,便于治疗过程的四手操作。

The rotary ceramic spittoon facilitates spitting by patients and facilitates four-handed operation.

组合脚开关

Multifunctional pedal switch


集合牙椅控制、吹屑气、手机水、冲盂、漱口多种功能一体,方便医生操作使用。

It integrates dental chair control, blowing air, hand piece water, flushing of the cuspidor and gargling water, convenient for doctors.



LED口腔灯

LED operation light

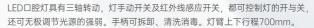

700mm

LED operation light has three-axis rotation, manual switch and infrared sensor switch, which can control the on-off and adjust the intensity. The handle can be disassembled, cleaned and disinfected. Vertical distance of the light arm is 700mm.



封闭电源

Enclosed power source

按欧洲安全规定要求,把高压部分封闭,防潮防泼溅,防止电磁干扰,提供稳定电源输出。

新引领未来

According to the requirements of European safety regulations, the high-voltage part is enclosed, moisture-proof and splash-proof, prevent electromagnetic interference and provide stable power output.


低压加热器

Low-voltage heater

低压加热器采用直流24V外加热方式 加热,避免了因在加热棒上附着水垢而改 变水质及加热效果,还可防干烧。

The low-voltage heater is heated by a DC 24V

external heating method, which avoids changing the water quality and heating effect due to the scale adhering to the heating rod and prevent dry burning as well.

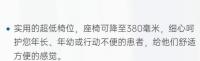

The practial super-low seat position 380mm, makes a convenient way for the elder, child and handicapped.

稳固、耐用的宽牙椅,靠背腰部凸型设计,给患者腰部良好支撑,特别是在仰卧时对整个背柱提供良好的承托。

The stable, durable and wide patient chair with lumbar support comforts the patient especially when he/she's slying on it and get good support on his/her back.



FOSHAN JOINCHAMP MEDICAL DEVICE CO., LTD.

地址:广东省佛山市禅城区科洋路1号(南庄罗格园高新工业区) 邮编:528000 总机:0757-82013008 传真:0757-82013009

销售部电话: 0757-82013003 售后服务部: 0757-82013012 出口部电话Export Dept: +86-757-82013002

WhatsApp: +86-13929980749 出口部邮箱: sales@chinadentalunit.com (Export Market)

中文网站: www.fszhongchuang.com 外文官网: www.zc-dentalunit.com

